Document made available under the Patent Cooperation Treaty (PCT)

International application number: PCT/IL05/000303

International filing date: 17 March 2005 (17.03.2005)

Document type: Certified copy of priority document

Document details: Country/Office: US

Number: 60/555,979


Filing date: 25 March 2004 (25.03.2004)

Date of receipt at the International Bureau: 09 June 2005 (09.06.2005)

Remark: Priority document submitted or transmitted to the International Bureau in

compliance with Rule 17.1(a) or (b)



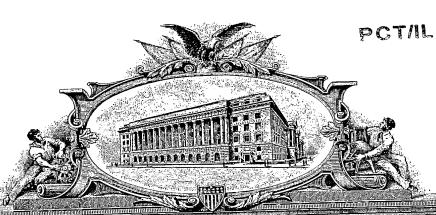

HIDE ONTO BED STRABES DEVINO BER COM

<u>TO AULTO WHOM THUSE; PRESENTS SHATU, COMER</u>

UNITED STATES DEPARTMENT OF COMMERCE

United States Patent and Trademark Office

May 03, 2005

THIS IS TO CERTIFY THAT ANNEXED HERETO IS A TRUE COPY FROM THE RECORDS OF THE UNITED STATES PATENT AND TRADEMARK OFFICE OF THOSE PAPERS OF THE BELOW IDENTIFIED PATENT APPLICATION THAT MET THE REQUIREMENTS TO BE GRANTED A FILING DATE UNDER 35 USC 111.

APPLICATION NUMBER: 60/555,979

FILING DATE: March 25, 2004

By Authority of the

COMMISSIONER OF PATENTS AND TRADEMARKS

P. R. GRANT

Certifying Officer

PTC/SB/16 (01-04)

Approved for use through 07/31/2006. OMB 0651-0032

U.S. Patent and Trademark Office; U.S. DEPARTMENT OF COMMERCE

Under the Paperwork Reduction Act of 1995, no persons are required to respond to a collection of information unless it displays a valid OMB control number.

PROVISIONAL APPLICATION FOR PATENT COVER SHEET

This is a request for filing a PROVISIONAL APPLICATION FOR PATENT under 37 CFR 1.53(c).

Express Mail Label No.

		IN	VENTOR(S)					
Given Name (first and middle [if any]) Family Name or Surname				(City a	Residence (City and either State or Foreign Country)				
ELAN		211	OM VI.		RAN	RAMAT GAN ISI		ISRA	EL
Additional inventors are be	sing named on the			perately numb	ered sheets a	ittached h	ereto		
	TIT	LE OF THE INVE							
A SOFT RIN		REVENTION		VAGWI					12
Direct all correspondence Customer Number:	to: CORF	RESPONDENCE AL	DDRESS	TR	LG AN	PRO	LAPS	e v	بز 2
OR	L		· · · · · · · · · · · · · · · · · · ·						s. F 97
Firm or Individual Name	ELAN	211	MD)						US. PTO 555979
Address	A HAIL	ANOT	ST						0/8/
Address									19
Gity	RAMAT	GAN	15	State		Zip	1526	48	
Country	ISRA		7	elephone		Fax 🚣	9+2-3	- 736	9268
		SED APPLICAT	ION PART	S (check all	that apply)				
		4			EVe\ Numbo				
Specification Number		<u></u>			D(s), Numbe				
Drawing(s) Number	of Sheets				ther (specify))			
,	eet. See 37 CFR 1.7								
METHOD OF PAYMENT	OF FILING FEES F	OR THIS PROVISION	ONAL APPL	CATION FOR	PATENT				
1 = 2	nall entity status. See	37 CFR 1.27.				FILING	G FEE		
A check or money	order is enclosed to	cover the filing fees	.				,	7	ļ
The Director is herby authorized to charge filing fees or credit any overpayment to Deposit Account Number.									
Payment by credit	card. Form PTO-203	38 is attached.				<u> </u>			
The invention was made United States Governme No. Yes, the name of the						ncy of the			
Respectfully submitted,	1		[Page 1 of 2		ate	21. :	3.04		
SIGNATURE	با ِ	(₀ 7		R	EGISTRATIC	NO.			•
TYPED or PRINTED NA	ME DR EV	9 N ZIV	MD	(it	appropriate) ocket Numbe				
TELEPHONE	972-53-	569928							
This collection of information to process) an application. gathering, preparing, and summent of time you require Trademark Office, U.S. Dep ADDRESS, SEND TO: Ma	USE ONLY is required by 37 CFR Confidentiality is gover to complete this form a partment of Commerce	FOR FILING A PR 1.51. The information ned by 35 U.S.C. 122 if application form to t and/or suggestions for P.O. Box 1450, Alex	ROVISIONAL n is required to 2 and 37 CFF the USPTO. r reducing this xandria, VA	o obtain or retain ? 1,14. This colle Time will vary de s burden, should 22313-1450. DO	a benefit by the ction is estimated pending upon be sent to the NOT SEND F	e public whated to take the individual Chief Informatics OR C	8 hours to out of case. Any mation Office OMPLETED	complete, y commer er, U.S. P. FORMS	including nts on the atent and TO THIS

If you need assistance in completing the form, call 1-800-PTO-9199 and select option 2.

. Solution the Paperwork Reduction Act of 1995, no person	ns are required to	respond to	U.S. Pate	nt and i	Tradomad Alle	mille ner	07/31/2006. OM	TOMBEDOE
EFF TO A NOW!	espond to a collection of information unless it displays a valid OMB control number. Complete if Known							
E FEE TRANSMIT	IAL	Applica	ation Nu	ımber				
for FY 2004		Filing Date		MARCH	MARCH 12 2004			
		First Named Inventor		r Dr Elan Z	Dr Elan Ziv			
Effective 10/01/2003. Patent fees are subject to ann		Examiner Name						
Applicant claims small entity status. See 37 CF	R 1.27	Art Unit						
TOTAL AMOUNT OF PAYMENT (\$) 80		Attorn	ey Dock	et No.				
METHOD OF PAYMENT (check all that ap	oply)			FEE C	CALCULATION	ON (conti	nued)	
Check Credit card Money Other	None 3.	ADDITIO	ONAL	FEES				
Deposit Account:		e Entity						1
Deposit	Fee		Fee Fo Code (ee \$)	Fee De	scription		Fee Paid
Account Number	105	51 130	2051	65 Su	urcharge - late fi	ling fee or o	ath	
Deposit Account	105	52 50	2052		urcharge - late p over sheet	rovisional fili	ing fee or	
Name The Director is authorized to: (check all that apply)	105	53 130	1053		on-English speci	ification		
Charge fee(s) indicated below Credit any over			1812 2,		or filing a reques	-		
Charge any additional fee(s) or any underpayment of fe	e(s) 180	04 920°	1804		equesting public xaminer action	ation of SIR	prior to	
Charge fee(s) indicated below, except for the filling fee	e 180	05 1,840*	1805 1		lequesting public examiner action	cation of SIR	after	
to the above-identified deposit account. FEE CALCULATION	125	51 110	2251	•	extension for rep	ly within first	t month	
1, BASIC FILING FEE	125	52 420	2252		extension for rep	_		
Large Entity Small Entity	129	53 950	2253	475 E	Extension for rep	within thin	d month	
Fee Fee Fee Fee Description Code (5) Code (5)	Fee Paid 12	54 1,480	2254	740 E	Extension for rep	oly within fou	rth month	
1001 770 2001 385 Utility filing fee	12	55 2,010	2255	1,005 E	Extension for rep	aly within fifth	month	
1002 340 2002 170 Design filing fee		01 330	2401		Notice of Appeal			
1003 530 2003 265 Plant filing fee 1004 770 2004 385 Reissue filing fee		02 330 03 290	2402		Filing a brief in s Request for oral !		appeal	
, , , , , , , , , , , , , , , , , , , ,		51 1,510			etition to institut	•	se amesedina	
SUBTOTAL (1) (\$) 80		52 110	2452		Petition to revive	•		
	14	53 1,330	2453	865 F	Petition to revive	- unintentio	nai	
2. EXTRA CLAIM FEES FOR UTILITY AND Fee from	10	01 1,330	2501	665 L	Utility issue fee (or reissue)		
Total Claims below		02 480	2502		Design issue fee	1		
Independent -3**= X		i03 640 I60 130	2503 1460		Plant issue fee	^		
Multiple Dependent		907 50	1807		Petitions to the (Processing fee t			
Large Entity Small Entity	8	180	1806		Submission of In			
Fee Fee Fee Fee Description Code (\$)	1	21 40	8021	AD F	Recording each	patent assign	nment per	
1202 18 2202 9 Claims in excess of 20	18	109 770	2809	, b	property (times n Filing a submissi			
1201 86 2201 43 Independent claims in ex	cesa of 3			((37 CFR 1.129(e	3))	-	
1203 290 2203 145 Multiple dependent claim 1204 86 2204 43 ** Reissue independent		310 770	2810	385 1	For each addition examined (37 Cl	nal invention FR 1.129(b)	to be	
over original patent		801 770	2801		Request for Co			
1205 18 2205 9 "Reissue claims in exce and over original pater	ess of 20 18 It	802 900	1802	900	Request for export of a design application	pedited exam	nination	
SUBTOTAL (2) (\$)		ther fee (s	pecify) _					
**or number previously paid, if greater, For Reissues,	see above	Reduced by	y Basic F	iling Fe	e Paid S	UBTOTAL	(3) (\$)	
SUBMITTED BY						(Complete (
Name (Print/Type) Dr Elan Ziv MD		Registra (Attorne)	ation No.	T			+9723736926	1
Signature . \` (33		CENSIN-			Date	March 14, 20	

WARNING: Information on this form may become public. Credit card information should not be included on this form. Provide credit card information and authorization on PTO-2038.

This collection of information is required by 37 CFR 1.17 and 1.27. The information is required to obtain or retain a benefit by the public which is to file (and by the USPTO to process) an application. Confidentiality is governed by 35 U.S.C. 122 and 37 CFR 1.14. This collection is estimated to take 12 minutes to complete, including gathering, preparing, and submitting the completed application form to the USPTO. Time will vary depending upon the individual case. Any comments on the amount of time you require to complete this form and/or suggestions for reducing this burden, should be sent to the Chief Information Officer, U.S. Patent and Trademark Office, U.S. Department of Commerce, P.O. Box 1450, Alexandria, VA 22313-1450. DO NOT SEND FEES OR COMPLETED FORMS TO THIS ADD SEND TO: Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

If you need assistance in completing the form, call 1-800-PTO-9199 and select option 2.

A SOFT RING FOR PREVENTION OF VAGINAL WALL NECROSIS IN TREATMENT OF PELVIC ORGAN PROLAPSE.

The present invention relates generally to the field of treatment of pelvic organ prolapse in female patients. More specifically, the present invention relates to a disposable device for use in the treatment of pelvic organ prolapse in women. The unique invention describes a soft ring covering a vaginal disposable device which is inserted and removed in a no-self-touch technique, by the patient herself, using a disposable applicator.

10

20

25

30

35

5

Inventor: Dr Elan Ziv, MD OBGYN, Urogynecologist

Background of the Invention

15 General Data

Definition

Pelvic organ prolapse is defined as a condition in which vaginal wall support is lost, and various pelvic organs prolapse into the vagina. This is a very bothering condition, though in most cases it is not a dangerous one. It might appear alone or in combination with urinary stress incontinence.

POP (Pelvic Organ Prolapse) is a very common condition in females, in which various organs, surrounding the vagina, prolapse into it. The reasons for such prolapse are numerous, mainly because the vagina has a very low tensile strength within its walls due to damage to muscles, nerves, fascias etc. There might also be a change within the collagen content, thereby causing a weaker pelvic floor. According to the older definition, prolapse may be divided into five categories, according to the organ that is sagging down (urethra, bladder, uterus, rectum and the pouch of Douglas (small bowl)), and to three grades according to the amount of descent (within the vagina, at the entrance to the vagina, protruding out of the vagina). There might be a combination of various organs prolapse at the same time, with different level of descent. The newer classification (POP-Q) takes into account other factors, such as location of the prolapse and the distance from the entrance of the vagina. There is very little data regarding the prevalence of the problem but it seems to be age related, and older women have much higher tendency towards developing some form of pelvic organ prolapse.

[&]quot;A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse"

Dr Elan Ziv, MD OBGYN, Urogynecologist Page 1 of 8

High risk group for Pelvic Organ Prolapse

There are very few statistical & Epidemiological studies regarding the causes for the development of pelvic organ prolapse and high risk groups. The main risk factors are:

- Advanced age
- Pregnancy
 - Delivery. During labor, there are many tears, stretchings and hematomas within the muscles, the ligaments and the nerves of the pelvic floor. Instrumental delivery might be a further factor.
 - Number of deliveries
- Race. White women tend to have more prolapse than black or oriental women.
 - Menopause
 - Obesity
 - Hard labor.

15 Prevalence

20

It was found that the prevalence of severe pelvic organ prolapse in women with jobs in labor/factory positions was 11.1%, 6.6% among those working at home, 3.0% in professional/ managerial jobs, 3.3% in service, and 0.8% among technical/sales/clerical workers. Laborers/factory workers had significantly more severe prolapse than the other job categories (P < .001). Table 1 summarizes some article on the subject

	Samuelsson et al 1999	Bland et al 1999	Progetto Menopausa 2000	Swift 2000 U.S.	
	Sweden	U.S.	Italy		
Country		241	21,449	497	
V	487	43-55	Around menopause 1	8->70	
Age range	20-59	45-35			
Race		89%	100%	47%	
White	100%	9%		52%	
Bluck	54%	Unspecified	85.6%	93%	
arcus	4%	28%	None	28%	
lysterectumy	Above, to, beyond hymen	ics	Utenis only	ICS	
Staging technique	Above, to, beyond hymen	100	Baden 0, 1, ≥2		
Stage		73%	94.5%	6.4%	
ō ₩	71.2%		3.6%	43.3%	
1	28.8%	23%	1.9%	47.7%	
2	2%	4790		2.6%	
3	0	0	_ ,		
4	0	0		Routine gynecolo	
Source of subjects	76% of every 3 year screen participants	Responded to ud for soybean supplement study	268 first-level menopause clinics in Italy.	gie examination	

TABLE 1

[&]quot;A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse"

Page 2 of 8

Page 2 of 8

Anatomy and pathophysiology

Understanding the female genital anatomy is extremely important in order to understand the basis of this invention.

- FIG. 1A is an anatomical frontal cut section of the vagina
- FIG. 1B is a front view of vagina creating the shape of the letter H. 5
 - FIG. 1C is a side view of the vaginal diameters.
 - FIG. 2 is a side view of the female pelvis.

10

15

.20

25

30

The vagina is a hollow organ with multiple shapes and diameters at different depths. The sidewalls (FIG.1A) have thin longitudinal and circular muscles (8), externally wrapped by the pelvic fascia (10), which tends to become weakened, stretched and attenuated over the years and allow for sagging and prolapse of different pelvic organs into and through the vagina. Internally the vagina is covered with the mucosa (6).

Although potentially being a tube, the vaginal cavity is almost always obliterated with the anterior wall (12) and the heavy bladder on it collapse and fall into the vaginal lumen. The posterior wall is pushed anteriorly by the bowel (14), thereby creating the shape of the letter H (FIG.1B). Vaginal diameters, either lateral or antero-posterior (FIG.1C, FIG.2) are never the same along vaginal axis. Lateral diameter is usually shorter closer to the entrance (18) gradually becoming longer internally (16). The same applies for the antero-posterior diameter, while diameter close to the entrance (44) is usually shorter than the one behind the perineal body (48). The vagina may therefore be looked at as a funneled shape organ.

There are certain defined organs in direct contact with the vagina or within short proximity. The bladder (36) is a hollow, sack like organ, in which the urine is accumulated prior to it's expulsion outside the body through the urethra (42). The bladder is located behind the pubic bone (38) and rests on the middle third of the vagina. The urethra is a short muscular pipe, 25-35 mms in length, in direct contact with the bladder at the bladder neck (40) resting on the lower third of the anterior vaginal wall, in direct contact with the pubic bone. The uterus (30) is a pear shaped organ that has a lower part, the cervix (32) which protrudes into the dome of the vagina. The body of uterus is an abdominal organ, which, in most cases bands forward on the bladder. The cervix, while protruding into the vagina creates 4 fornices - 2 laterals which are of the same dimensions, one smaller anterior (34), and one posterior (22), which is actually

[&]quot;A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse Page 3 of 8 Dr Elan Ziv, MD OBGYN, Urogynecologist

behind the cervix. The posterior fornix has direct contact with the pouch of Douglas (20) which is the deepest part of the abdominal cavity, on top of the vagina. The posterior surface of the vagina is in contact with the rectum (24), until the area where a thick muscular body, the perineal body (26) which is part of the pelvic floor and the sphincter of the anus (28).

Prolapse of adjacent organs into the vagina is quite common, in variable degrees. The reasons for such a prolapse are mainly damage to the endopelvic fascia which surrounds the organs and keeps them in the right position, pelvic floor muscular damage, and neural damage. Prolapse of the urethra, named urethrocel, is one of the several possibilities for pelvic organ prolapse around the vagina. The most frequent prolapse is of the bladder, named cystocele, as a result of its weight. Prolapse of posterior fornix, with bowl inside is called enterocele, and prolapse of the posterior vaginal wall with the rectum in it is termed rectocele. When the ligaments that hold the uterus in place weaken, a uterine descent occurs. In cases where the uterus has been removed, and the vagina is dome shaped, a vault prolapse may occur.

Current treatment

5

10

15

.20

30

At present - there are 2 ways of dealing with Pelvic Organ Prolapses:

- 1. Surgical vaginal or abdominal
- 2. Use of vaginal devices (pessaries) that are inserted into the vagina and mechanically reduce the prolapse by pushing the wall aside and upwards.

Efforts to avoid surgical procedures have resulted in the development of a number of non-surgical vaginal devices, inserted into the vagina by the surgeon or the patient.

- Vaginal devices are well known for their tremendous diversity in shapes and sizes.

 These devices are meant to prevent prolapse of the vaginal walls, with different location of pressure application.
 - Some of these devices tend to block all flow of urine from the bladder. Therefore, when a patient needs to urinate, the device must be removed from the vagina or must be collapsed to remove the pressure applied against the bladder neck. Trying to solve this problem, vaginal devices were developed in special shapes, without completely blocking the bladder neck so that the patient may urinate with the device in place (for example, US patents 6189535 B1, 6158435, 5894842, 5611768). These devices,

[&]quot;A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse"

Dr Elan Ziv, MD OBGYN, Urogynecologist Page 4 of 8

however, are generally large and intrusive and, therefore, are uncomfortable to wear, with low patient satisfaction and compliance. They are also relatively expensive, and therefore designed to be reusable.

5 There are several drawbacks of existing vaginal devices:

10

15

20

30

- 1. Most of them are intended to be reusable, hence they are made as resilient large bodies, made of plastic, hard rubber, or other such materials, in order to preserve their shape and function for a long time.
- 2. Insertion of large noncompliant bodies is sometimes difficult, painful or unpleasant, sometimes necessitating a medical practitioner. Removal causes same unpleasantness.
- 3. Most reusable devices are meant to remain in the vagina for prolonged periods of time, thereby causing irritations, pressure ulcers, infections, foul smelling discharge, etc.
- 4. Such reusable devices have a long standing bad reputation among patients and medical practitioners for being unpleasant, causing infectious discharge and foul odor, and being associated with disability and old age.
 - 5. Some reusable devices are meant to be inserted daily by the patient, and to be removed after several hours, by means of pulling a string or with a finger, to be cleaned for re-use, and to be kept in certain conditions prior to following insertion. Some patients are reluctant to touch their selves in such intimate parts of their bodies, or disgusted to clean the device, hence their reluctance to use it.

Vaginal devices may cause two main side effects:

- 25 Infections-any foreign material, anywhere in the body, may become infected by several kinds of organisms, and cause formation of a foul smelly discharge. Reusable devices are certainly prone to cause such infections, but also, less frequently, disposable ones. In order to prevent such infections, devices should be:
 - disposable
 - used for a limited length of time, in order to prevent organism from growing, and to prevent damage, by vascular pressure, to the vaginal wall.

[&]quot;A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse"

Dr Elan Ziv, MD OBGYN, Urogynecologist

Page 5 of 8

- made in a way that will allow vaginal or cervical normal secretions to flow out of the body
- made of certain known materials that would not permit growth of organisms.
- Pressure Necrosis prolonged pressure on the vaginal sidewalls may cause pressure on blood vessels with resultant necrosis, bleeding and infections. This might be prevented by using disposable devices only, for a predetermined length of time.

10

15

[&]quot;A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse"

Dr Elan Ziv, MD OBGYN, Urogynecologist Page 6 of 8

The invention

10

15

20

25

30

The invention will now be described with reference to accompanying drawings:

FIG. 2 is a side view of the female pelvis.

FIG. 3A is a top view of the invention partially covering a vaginal device.

5 FIG. 3B is a side view of the invention partially covering a vaginal device.

FIG.4A is a sectional view of the invention. (Section A-A)

FIG.4B is a sectional view of the invention. (Section B-B)

The device (FIG. 3) is intended to be applied into the vagina with an applicator much like the insertion of a regular menstrual tampon. After insertion, the device should expand significantly to the predefined shape and size, thereby exerting predefined appropriate pressure on both lateral vaginal walls, pushing them aside. The apex of the vagina shall be pushed upwards at the same time. That expansion of the device shall eventually create linear stretching of the anterior & posterior vaginal walls, while creating a new shape of intra-vaginal hollow (rectangle). The device may be left inside the vagina for several hours. Removal for disposal will occur while pulling a string and collapsing the device to a much smaller size.

The vaginal device (FIG. 3A+B), looks like a round ring (50+52) and a string (54). It is made of elastic materials, such as silicon, polyurethane, etc. This flexible material gives the device the ability to be folded into an applicator, and to be pushed out easily within small dimensions into the vagina.

On its outer surface, the device (50) is partially covered by an external soft layer (52) made of a much softer material, such as sponge rubber, etc., which makes it much more comfortable to wear. The sponge is also reducing the pressure on the vagina's sidewalls and helps preventing pressure necrosis. This prevention of pressure takes place due to the spongy texture of the device, thus blood flow within the vaginal walls is not halted, and pressure necrosis is less likely to happen.

The presence of such a spongy layer further reinforce the friction between the device and the walls of the vagina, to further prevent unwanted movements or expulsion of the device.

The spongy layer is situated on the external aspect of the device (FIG 3A+3B). In fact, only part of the external wall of the device (58) is covered (56) as is shown in FIG 4.

[&]quot;A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse"

Dr Elan Ziv, MD OBGYN, Urogynecologist Page 7 of 8


Alternative embodiments of the invention.


5


- The external soft layer may be made of several soft materials
- Degree of coverage of the device with the spongy layer might be altered according to shape, pressure and size
- Thickness of the spongy layer might also be changed.

"A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse"


Dr Elan Ziv, MD OBGYN, Urogynecologist Page 8 of 8

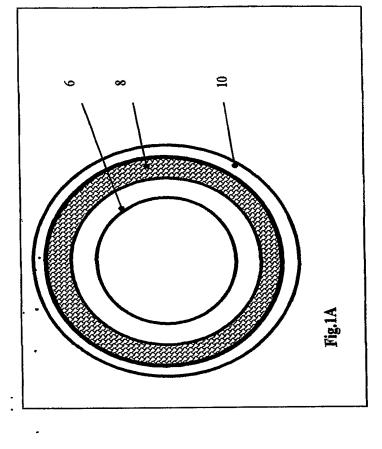


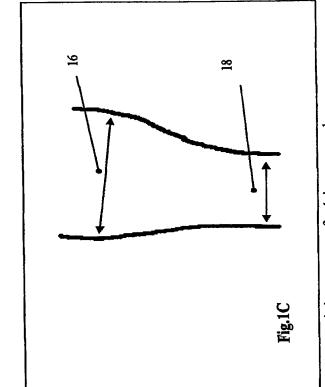

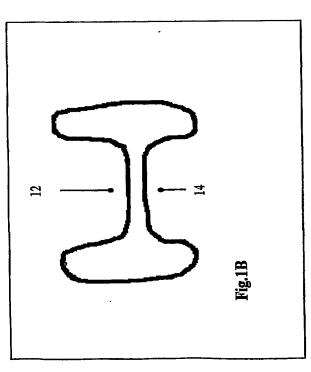

A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse Dr Elan Ziv, MD OBGYN, Urogynecologist

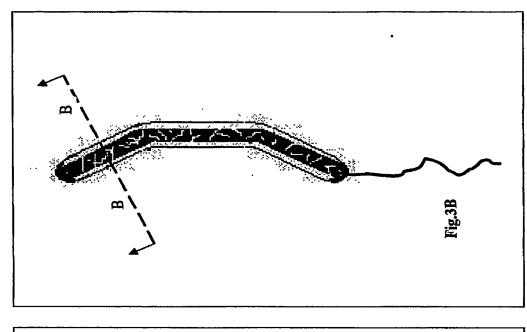

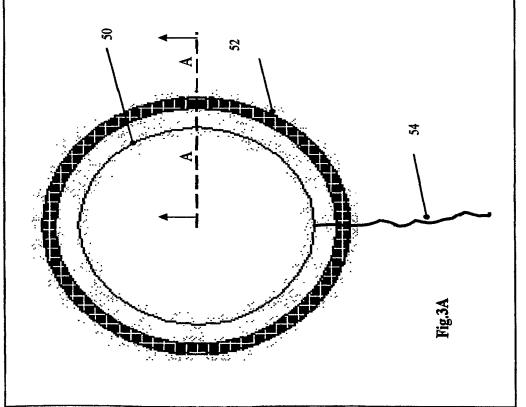

A soft ring for prevention of vaginal wall necrosis in treatment of pelvic organ prolapse Dr Elan Ziv, MD OBGYN, Urogynecologist

